



FOR FLUID MILK SAMPLES With 2% FAT

Catalogue No.: KTSP-71051 96 tests

Enzyme immunoassay kit for the quantitative determination of Vitamin-A in fluid milk containing 2% milk fat

For *in vitro* quantification use only.

I. PROPRIETARY NAME

The *VitaKit A™* from SciMed Technologies Inc., Catalogue No. KTSP-71051 contains sufficient material to assay 96 tests.

II. APPLICATION AND INTENDED USE

Dairy milk is fortified with vitamins A & D, as milk has become the major source of these vitamins for human beings. Regulatory agencies have set standards specifying the amount of vitamins A and D to be added to milk products. The methodology has been designed to extract Vitamin A (all trans retinol) from milk fat, and to directly quantify the Vitamin A content using an ELISA based assay. Other methods that can detect vitamins in dairy milk are time consuming and require specialized laboratory equipment and trained personnel. The *VitaKit A™* provides materials for the quantitative measurement of vitamin-A in fluid milk containing 2% milk fat for routine quality control purpose. This assay is intended for *in vitro* quantification only.

III. PRINCIPLE OF THE METHOD

The vitamin A ELISA test is based on the principle of a sandwich enzyme-linked immunosorbent assay. The assay system utilises immobilization of a monoclonal antibody directed against a distinct antigenic determinant on the vitamin-A molecule (retinol) on a solid phase. A second Anti-vitamin-A monoclonal antibody conjugated to horseradish peroxidase (HRP) is used as the detecting antibody in the assay. The test sample is allowed to react sequentially with the two antibodies, resulting in the vitamin-A molecules being sandwiched between the solid phase and enzyme-linked antibody. After the incubation, the wells are washed with distilled water to remove all

unbound labelled antibodies. TMB solution is added as a substrate, resulting in the development of a blue color. The reaction is stopped with the addition of 0.2M H₂SO₄, changing the color to yellow. The concentration of vitamin A is directly proportional to the color intensity of the test sample. Absorbance is measured spectrophotometrically using a plate reader at 450 nm.

IV. REAGENTS SUPPLIED WITH KIT

Storage : 2 °C – 8 °C
Stability : refer to expiration date on reagent labels

- [SORB]** Anti-Vitamin-A antibody coated wells : **[REF]** CW-71051: 96 wells with Anti-Vitamin-A Mab immobilized on the wells.
- [CONJ]** **[ENZ]** Anti-Vitamin-A antibody conjugate with HRP **[REF]** EC-71051: one (1) vial containing 0.1 mL of concentrated Anti-Vitamin-A conjugated with HRP, in a stabilizer solution.
- [CAL]** **[1-5]** Vitamin-A Standard **[REF]** WSC-71051: Standards prepared in hexane: 0.0, 0.055, 0.11, 0.22, 0.45 IU/mL. Content is 0.5 mL per vial.
- [CONTROL]** **[1]** Control 0.3 IU/mL **[REF]** QC-71051: 0.5 mL per vial.
- [CONTROL]** **[2]** Control 0.1 IU/mL **[REF]** QC-71052: 0.5 mL per vial.
- [BUF]** Reaction Buffer **[REF]** RB-71051 one (1) vial containing 13 mL of protein based (milk) buffer with thimerosal as preservative.
- [SUBS]** **[TMB]** Enzyme substrate **[REF]** ES-71051: one (1) vial containing 11 mL of TMB solution.
- [CONJ]** **[DIL]** Conjugate Diluent **[REF]** CD-71051: one (1) vial containing 11 mL of glycine based Buffer.
- [H₂SO₄]** Stopping solution **[REF]** SS-71051: one (1) vial containing 12 mL of 0.2 M sulfuric acid.
- [STAB]** Stabilizer **[REF]** S-71051: one (1) vial containing 0.15 mL of stabilizer.

V. EQUIPMENT & MATERIAL REQUIRED BUT NOT PROVIDED

- ✓ Precision pipettes with disposable tips
- ✓ 8 channel pipette (100-200 µL) with disposable tips
- ✓ Plate shaker set at 180 ± 10 rpm
- ✓ Microplate reader with filter at 450 nm
- ✓ Microplate washer
- ✓ Centrifuge
- ✓ Vortex Mixer
- ✓ Deionized or distilled water
- ✓ Absorbent paper
- ✓ 95% ethanol (denatured)
- ✓ Potassium hydroxide (KOH) pellets
- ✓ Hexane
- ✓ 10 mL screw capped glass tubes
- ✓ 1 or 2 mL screw capped amber glass vials

VI. REAGENT PREPARATION

- Bring all reagents to room temperature before use (22 °C ± 2 °C), except enzyme conjugate concentrate **[CONJ]** **[ENZ]** (EC-71051) which should be at kept 2 °C - 8°C
- Enzyme conjugate concentrate (EC-71051) should be diluted as indicated on the bottle (label) with conjugate diluent (CD-71051) according to the number of wells used.** Mix the enzyme conjugate concentrate by pipetting 2-3 times with a pipette tip before diluting with diluent. Add required amount of enzyme conjugate concentrate to the conjugate diluent and mix thoroughly before use. Diluted conjugate can not be stored and should be prepared fresh for each assay run.

Handling notes:

Do not mix materials from different kit lots. Bring all reagents to room temperature before using except Anti-Vitamin-A antibody conjugate with HRP. Use a clean disposable pipette tip for addition of each different sample and reagent to avoid cross-contamination. Use only glass vials for the extraction of vitamins. Perform extraction steps in amber light or clear glass vials should be wrapped in aluminium foil if working under fluorescent light. Use a stop watch as timing is critical.

Prepare a standard curve for each run. Do not use data from previous runs. Cap all Vitamin-A calibrators and Vitamin-A extracted specimens immediately after loading onto ELISA plate. This will allow the reference calibrators and extracts to be used more than once if desired. Load all extracted specimens and reference calibrators within 5 minutes and accurately onto the ELISA strips to limit variations in evaporation time between the first and last well loaded. All hexane steps should be performed under the fume hood. Avoid contact with eyes and skin.

VII. EXTRACTION PROCEDURE

Fluid milk samples should be stored at 2°C -8°C. Bring fluid milk container to room temperature. Rotate slowly at least 10 times without foaming. Pipette 1 mL of milk sample using a clean pipette into the corresponding glass tube. Store the remaining milk sample at 2°C - 8°C.

A. Milk with 2%M.F.

- Label 10 mL screw capped glass tubes and pipette 1 mL of milk in corresponding tube. Perform extraction steps in amber light or clear glass vials should be wrapped in aluminium foil if working under fluorescent light.
- Pipette 15 µL of ethanol into each tube, Do not cap the tubes. Shake for 20 seconds in the dark.
- Add 0.60 g of KOH into milk with 2% milk fat., Gently mix for 2 minutes in the dark. Cap the tubes and incubate at room temperature for 4 minutes in the dark. Shake vigorously for 2 minutes. Repeat 4 minutes incubation and 2 minutes vigorous shaking 2 more times (totals 12 minute incubation and 6 minute shaking).
- Pipette 2 mL of hexane into above solution. Cap and shake vigorously for another 2 minutes in the dark.
- Centrifuge each tube at room temperature for 5 minutes at 2400 rcf. After centrifugation, handle tubes carefully. The upper organic phase must be perfectly clear and well separated.
- Label 1 or 2 mL screw capped amber coloured glass vials and add 300 µL of hexane into each vial for dilution of vitamin-A extract. Transfer 100 µL of vitamin-A extract in corresponding

amber coloured glass vials (100 µL of organic phase + 300 µL of hexane). The amber coloured glass vials, which contain 8-fold diluted vitamin A extract, must be capped well and should be assayed immediately. If necessary dilute the extract with more hexane.

Table I

	Steps	2%M.F.	Conditions
Saponification and extraction	Fluid Milk	1 mL	Run water control only for skim milk.
	Ethanol 95%	15 µL	Gently mix for 20 seconds.
	KOH (g)	0.6	Gently mix for 2 minutes in the dark.
Extract dil. with hexane	Hexane	2 mL	Shake vigorously for 2 minutes in the dark and centrifuge at 2400 rcf for 5 minutes.
	Upper organic phase	100 µL	The diluted extract in screw capped amber coloured glass vial should be assayed immediately.
Hexane	300 µL		

VIII. ASSAY PROCEDURE

Refer to the assay procedure, Table II. Standards, specimens and controls should be assayed in duplicate. Secure the desired number of coated wells **[SORB]** in the holder.

- Pipette 10 µL of calibrators **[CAL]** **[1-5]**, controls **[CONTROL]** **[1]**, **[CONTROL]** **[2]** and diluted extracted specimens into the corresponding wells.
- Shake the wells 6 minutes on a plate shaker (180 ± 10 rpm) at room temperature (22 ± 2°C) to evaporate hexane
- Pipette 120 µL of Reaction Buffer **[BUF]** into each well. Mix gently for 30 seconds. Place opaque lid or adhesive cover over the strips.
- Incubate for 30 minutes in the dark on the plate shaker (180 ± 10 rpm) at room temperature (22 ± 2°C).
- Wash four times with distilled water using Microplate washer. Manual washing may also be used: with wash bottle or using multi-channel pipette, add 300 µL of distilled water in each well in each wash cycle. Care should be taken to avoid spillage of distilled water into adjacent wells. After the wash, completely decant the water by tapping the plate against absorbing paper until no trace of water is visible on the paper.
- Mix the freshly diluted Anti-Vitamin-A conjugate-HRP. Pipette 100 µL of freshly diluted Anti-Vitamin-A conjugate-HRP **[CONJ]** **[ENZ]** in each well. Mix gently for 20 seconds. Place opaque lid or adhesive cover over the strips.
- Incubate for 20 minutes in the dark at room temperature (22 ± 2°C).
- Wash 6 times (refer to step no. 5).

9. Pipette 100 µL of TMB **SUBS** **TMB** (Substrate) into each well. Gently mix for 10 seconds.
10. Incubate 10 minutes in the dark at room temperature (22 ± 2°C).
11. Add 100 µL of the stopping solution **H₂SO₄**. Gently mix for 10 seconds.
12. Measure the absorbance at 450 nm using a microplate reader.

NOTE: READ THE ABSORBANCES IMMEDIATELY AFTER COMPLETING THE ASSAY.

TABLE II


Wells	Identification	Assay Volumer	Evaporate	Reaction Buffer		Dil. Conjugate		Substrate	Stop. Sol.	
A ₁ ,A ₂ B ₁ ,B ₂ C ₁ ,C ₂ D ₁ ,D ₂ E ₁ ,E ₂ F ₁ ,F ₂ G ₁ ,G ₂ H ₁ ,H ₂ etc...	0 IU/mL 0.055 IU/mL 0.11 IU/mL 0.22 IU/mL 0.45 IU/mL Control 1 Control 2 Sample extract etc...	100 µL		120 µL INCUBATE	DECANT & WASH	100 µL INCUBATE	WASH	100 µL INCUBATE	100 µL	READ AT 450 nm

IX. CALCULATIONS

The standard curve is used to determine the amount of Vitamin A in unknown sample. The standard curve is generated by plotting the average of O.D. (450 nm) obtained for each of the standard concentrations on the vertical (y) axis versus the corresponding standard concentrations on the horizontal (X) axis.

- Examine data for acceptance criteria with quality control guidelines.
- For unknown samples, multiply the values by a factor of 8 as per the stipulated assay conditions.
 - 3.33 I.U. of Vitamin A = 1 µg

EXAMPLE OF Vitamin-A STANDARD CURVE

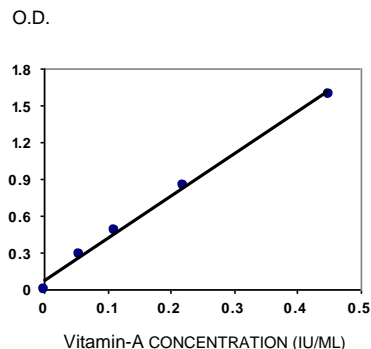

X. PERFORMANCE CHARACTERISTICS

Sensitivity: The range for this assay under the specified conditions is from 0.055 I.U./mL to 0.45 I.U./mL.

Limit of Detection: The minimum detection limit for this assay is 33 IU/100 mL or 0.33 IU/mL.

Precision & reproducibility: The relative standard deviation for interassay and intrassay was determined to be 8.6 % and 1.9% respectively.

Cross Reactivity and Specificity: The kit did not exhibit any significant cross reactivity with cholesterol and vitamin D3 and is specific for Vitamin A.

Standard curve Linearity Linearity was determined to be 0.996 (Average of six independent assays) with %RSD of 1.5%.

XI. Limitations of the procedure:

- a) Reliable and reproducible results will be obtained when the assay procedure is carried out with strict adherence to the procedure described within this package insert and good laboratory practice.
- b) A maximal total pipetting time of 5 minutes for calibrators, controls and specimens is suggested.

XII. QUALITY CONTROL

Good laboratory practice requires that quality control specimens be run with each calibration curve to check the assay performance.

XIII. SAFETY MEASURES

1. All materials in this kit should be used only for *in vitro* quantitative tests not involving internal or external administration of the material to humans or animals.
2. Respect laboratory quality controls rules.

3. Reagents are matched in each kit, and, therefore, reagents from different lot numbers should not be mixed.
4. This kit should not be used after the expiration date.
5. Optimal results will be obtained by strict adherence to this protocol.
6. The stopping solution contains sulfuric acid. This solution should be handled with caution, avoiding skin contact.
7. KOH pellets should be handled with caution, avoiding skin contact.
8. All hexane steps should be performed under the hood. Avoid contact with eyes and skin. Do not inhale.
9. Prior to assay, Bring all reagents to ambient temperature by allowing them to stand at room temperature (22 ± 2°C). Gently mix all reagents.

XIV. Technical Assistance

SciMed Technologies Inc. Suite #209, Advance Technology center 9650 – 20 Avenue Edmonton, Alberta Canada T6N 1G1
 Phone: 780-702-1509
 Fax: 780-702-0303
 scimed@scimedtechnologies.com
 www.scimedtechnologies.com

XV. LIST OF REAGENTS SUPPLIED WITH KIT

Cat. #	Description	KTSP-71051
CW-71051	Anti-Vitamin-A coated wells	96 wells
EC-71051	Anti-Vitamin-A conjugate HRP	1
WSC-71051	CAL 1- CAL 5	1
QC-71051	Control 1	1
QC-71052	Control 2	1
CD-71051	Conjugate Diluent	1
RB-71051	Reaction Buffer	1
ES-71051	Enzyme Substrate (TMB)	1
SS-71051	Stopping Solution	1
S-71051	Stabilizer	1
	English Protocol	1

XVI. SYMBOLS ON REAGENTS LABEL



SciMed Technologies Inc. does not offer a warranty of any kind, neither express nor implied, other than that this product meets all the quality control standards set by the SciMed Technologies and will, at its discretion, replace any product which proves to be defective in terms of workmanship or material. SciMed Technologies Inc. is not liable in any way for loss of profits or any damage, direct or otherwise, to persons or property, arising from the use of this product or these procedures.

